

pfmmedical

Quality and Experience

ASEPT[®]

- › Pleural & Peritoneal Drainage Systems
- › Drainage Bottles, Bags & Accessories
- › Insertion Kits


Offering patients the benefit of both gravity & vacuum bottle drainage systems for hospitals and out-patient care.

www.pfmmedicalusa.com



ASEPT Catheter and Insertion Kits

The ASEPT Pleural and Peritoneal Drainage Systems are designed for long-term, intermittent drainage of symptomatic, recurrent fluid accumulation.

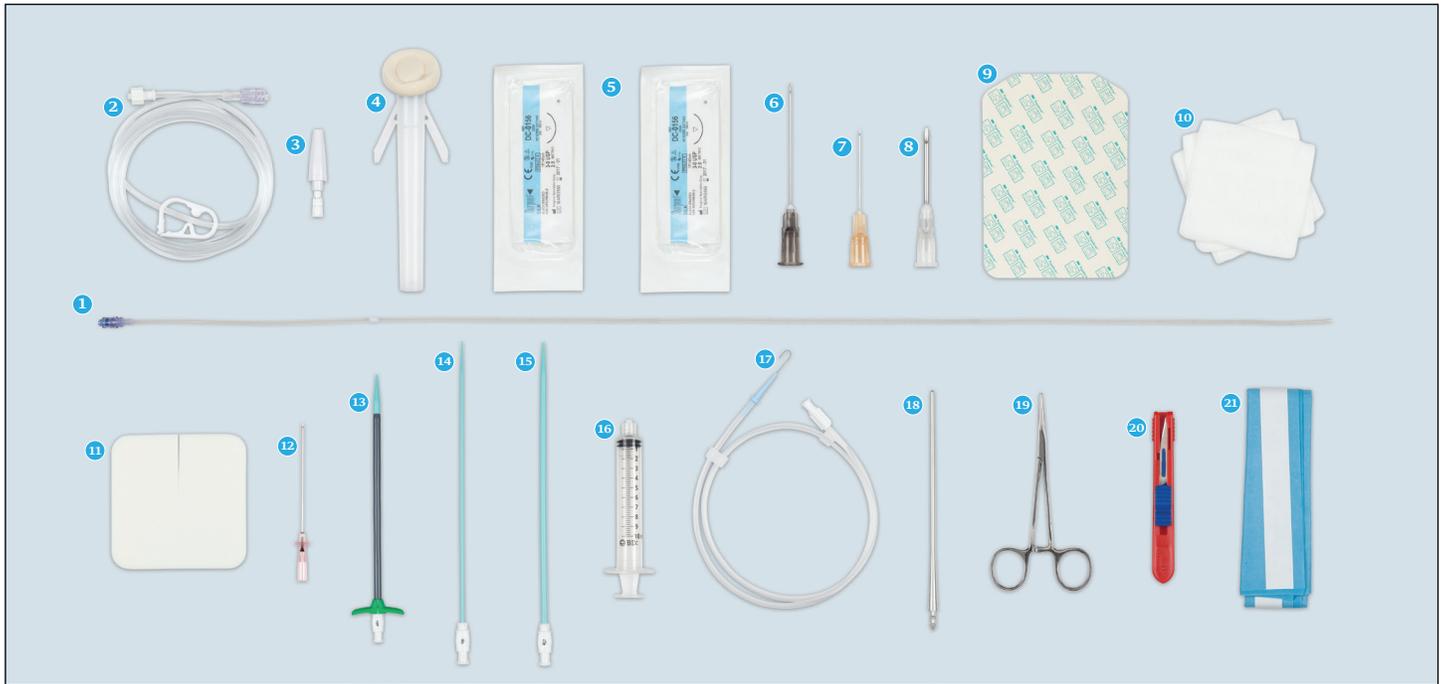

- | | | | |
|---|------------------------------------|------------------------|----------------------|
| 1 ASEPT Catheter | 6 16G x 1.0" (25.4mm) Needle | 12 Introducer Needle | 18 Tunneler |
| 2 ASEPT Drainage Line | 7 25G x 1.0" (25.4mm) Needle | 13 Tearaway Introducer | 19 Forceps |
| 3 5-in-1 Drain Line Adapter | 8 22G x 1.5" (38.1mm) Needle | 14 8F (2.67mm) Dilator | 20 Safety Scalpel |
| 4 Chloraprep® 10.5mL Vial | 9 Tegaderm® Self Adhesive Dressing | 15 12F (4.0mm) Dilator | 21 Fenestrated Drape |
| 5 Synthetic Suture, 2-0 Straight and Synthetic Suture, 3-0 Curved | 10 Gauze Pads, 4"x4" | 16 10mL Syringe | |
| | 11 Foam Catheter Pad | 17 J-Tip Guidewire | |


The primary component of the ASEPT Drainage system is the catheter, which has an aseptic, self-sealing valve attached that provides a secure and reliable connection. The valve can be connected to an ASEPT vacuum drainage bottle, drainage bag, or a stand-alone drainage line to connect to other fluid collection equipment.

Features:

- ▶ Secure connection
- ▶ Cap-less and Aseptic design
- ▶ Needle-free and swabable
- ▶ Easy to replace
- ▶ Self-sealing


ASEPT Drainage Bottles and Bags

The ASEPT Drainage kit configuration includes a 600mL or 1000mL vacuum bottle set, a 2000mL gravity drainage bag, and the procedure pack with all of the dressing change items needed to complete an at-home drainage procedure. A standalone drain line is available to connect to the ASEPT drainage catheter for drainage with other collection equipment during in-patient procedures.

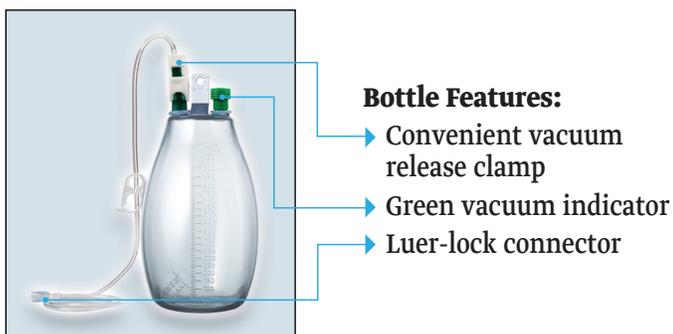

Bottle Features:

- ▶ Convenient vacuum release clamp
- ▶ Green vacuum indicator
- ▶ Luer-lock connector

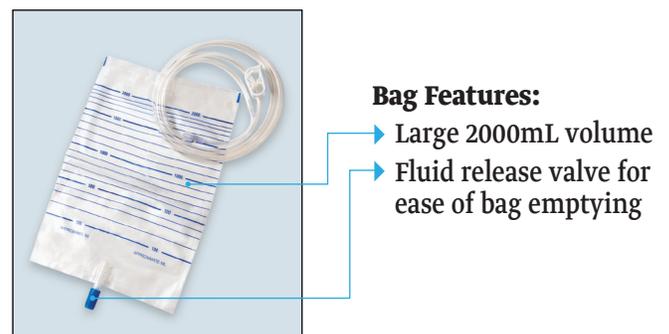

Bag Features:

- ▶ Large 2000mL volume
- ▶ Fluid release valve for ease of bag emptying



◀ **ASEPT Starter Kit**

The ASEPT Starter Kit is available to ease the transition for the patient from the hospital to their home. The ASEPT Starter Kit features an easy-to-carry box with 5 ASEPT Drainage bottles and Procedure Packs. The ASEPT Starter Kit provides peace-of-mind for the patient as the patient adapts to this new method of managing their own care at home.


ASEPT Education

The transition from the hospital to the home can sometimes be difficult for an ill patient. In an effort to ease some of the difficulties associated with this transition, pfm medical, inc. provides a plethora of educational materials to assure the patient that they are in good hands and perfectly capable of managing their care from the comfort of their own home.



◀ **Sales Support**

The first line of defense for any questions or concerns a patient may have about their new ASEPT drainage system will be the friendly and knowledgeable sales team at pfm medical, inc. All sales staff are thoroughly trained on the ASEPT system and equipped to assist any patient or caregiver with the knowledge and tools essential to perform the drainage procedure safely and comfortably.

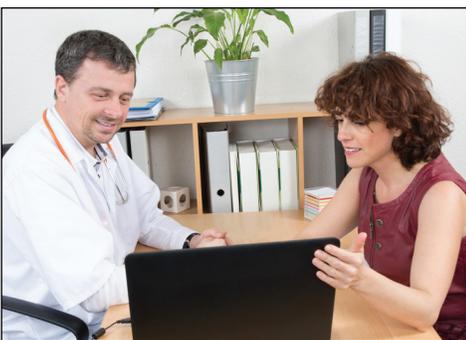

◀ **Training DVD**

Each ASEPT patient will receive a short, 9-minute video detailing the ASEPT drainage procedure. The video outlines each step that needs to be taken by the patient and/or caregiver to assure that the drainage procedure is completed accurately.

ASEPT Poster ▶

Each ASEPT patient will also receive an ASEPT Step-By-Step Guide. Like the ASEPT video, this guide displays the drainage procedure step-by-step. Designed for continual reference during repeated drainage procedures, the ASEPT Step-By-Step Guide reinforces what was taught in the ASEPT video and provides the patient or caregiver with a convenient reference to make sure no steps are missed.



ASEPT Accessories

pfm medical, inc. is proud to offer a full line of accessories ensuring that the physician or patient will have everything he or she needs to complete the drainage procedure. Our accessory portfolio contains the following items:

ASEPT Replacement Valve



ASEPT Drainage Line Set



ASEPT Procedure Pack



Ordering Information

ASEPT Insertion Kits

Item No.	Description	Box Qty
M7001	ASEPT Pleural Drainage System (Includes ASEPT catheter and insertion kit)	1 each
M7008	ASEPT Peritoneal Drainage System (Includes ASEPT catheter and insertion kit)	1 each
M7009	ASEPT Supplemental Insertion Kit (Includes 7F dilators & 16F valved peelable introducers)	10 each

ASEPT Drainage Bottles and Bags

Item No.	Description	Box Qty
M7002	ASEPT 600mL Drainage Kit (Includes vacuum bottle and Procedure Pack)	10 each
M7000	ASEPT 1000mL Drainage Kit (Includes vacuum bottle and Procedure Pack)	10 each
M7006	ASEPT 600mL Select Drainage Kit (Includes low vacuum bottle and Procedure Pack)	10 each
M7007	ASEPT 1000mL Select Drainage Kit (Includes low vacuum bottle and Procedure Pack)	10 each
M7060	ASEPT 2000mL Drainage Starter Kit (Includes gravity bag with drainage line - for Peritoneal drainage only)	10 each
M7058	ASEPT 2000mL Drainage Bag (Includes gravity bag with drainage line - for Peritoneal drainage only)	10 each
M7080	ASEPT 600mL Drainage Starter Kit (Includes vacuum bottle and Procedure Pack in patient box)	5 each
M7070	ASEPT 1000mL Bottle (Includes vacuum bottle - does not include Procedure Pack)	10 each
M7072	ASEPT 600mL Bottle (Includes vacuum bottle - does not include Procedure Pack)	10 each
M7081	ASEPT 1000mL Drainage Starter Kit (Includes vacuum bottle and Procedure Pack in patient box)	5 each

ASEPT Accessories

Item No.	Description	Box Qty
M7040	ASEPT Procedure Pack (Includes all necessary dressing change items)	10 each
M7004	ASEPT Replacement Valve	5 each
M7005	ASEPT Drainage Line Set (Includes drainage line set and 5-in-1 adapter)	10 each

Chloraprep is a registered trademark of Carefusion 2200, Inc. Tegaderm is a registered trademark of 3M.

www.pfmmedicalusa.com



pfm medical, inc.
 1916 Palomar Oaks Way, Suite 150 • Carlsbad, CA 92008
 Ph: 760.758.8749 • Fx: 760.758.1167
 customerservice@pfmmedicalusa.com • www.pfmmedicalusa.com

©2018 PFM. All rights reserved.
 Unless otherwise noted, PFM Medical, the PFM logo, the color GREEN of the three-dimensional configuration of the accordion vacuum indicator, and ASEPT are registered trademarks of pfm medical, inc.