

pfmmedical

Quality and Experience

JET-VAC®

- › Procedure Trays
- › Centesis Catheters
- › Drainage Bottles

Efficient, safe and easy-to-use short term percutaneous drainage systems for hospitals.

www.pfmmedicalusa.com


**Safety Tray
Now
Available!**



Jet-Vac Centesis Catheters, Jet-Vac Procedure Trays and Jet-Vac Bottles provide a complete drainage system for physicians, nurses, and hospitals. The Jet-Vac drainage bottle is a disposable vacuum device used to aspirate, remove, or sample body fluids.

Jet-Vac Centesis Catheters

- ▶ **Large circular side holes designed to increase fluid drainage**
- ▶ **Tapered catheter tip**
 - Provides a smooth transition
 - Helps minimize patient discomfort
- ▶ **HDPE material**
 - Reduces catheter collapsing
 - Promotes a smooth catheter insertion
- ▶ **Clear needle hub with catheter lock**
 - Increases fluid visibility
 - Ergonomic design provides a secure grip
 - Prevents needle/catheter separation
- ▶ **Plastic vent plug**
 - Allows for air to exit
 - Provides added safety for physician
- ▶ **Echogenic needle**
 - Ensures needle visibility during ultrasound placement

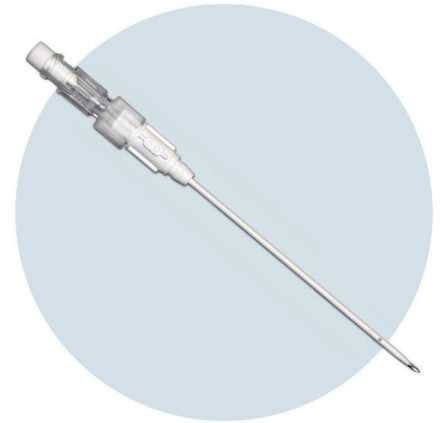

Jet-Vac Procedure Trays

- ▶ **All your drainage needs in one location**
 - Allows for drainage procedure to be performed anywhere with minimal supplies

Tray includes:

- | | |
|--|-------------------------------|
| 1 - Fenestrated Drape | 6 - Gauze Pads (4 x4) |
| 1 - ChloroPrep® 3 mL with Tint | 1 - #11 Safety Scalpel |
| 2 - Lidocaine HCL 1% 5 mL Ampules | 1 - 60 mL Syringe |
| 2 - Lidocaine Labels | 2 - 10 mL Syringe |
| 1 - Bandage Dressing | 1 - CSR Wrap |
| 2 - Specimen Vial, Clear 10 mL | 1 - Procedure Tray |
| 1 - Large Bore Stopcock | 1 - 5-in-1 Adapter |
| 2 - Filter Straws 1.75" Hydrophilic | 1 - Spike |
| 1 - 25G, 1.5" Needle* | 1 - Drainage Line 48" |
| 1 - 22G, 1.5" Needle* | 1 - Jet-Vac Centesis Catheter |
| 1 - Needle Holder (not included in Safety Kit) | |



*Safety Needle option available, see ordering information

Jet-Vac Bottles

- ▶ **Needle-Free Luer Locking Connection**
 - Removes need for spikes or needles to access the bottle
 - Sample fluids aseptically
- ▶ **Shatter-Proof Plastic Design**
 - Reduces risk for hospital personnel and patients by avoiding glass breaks
 - Reduces hazardous waste disposal costs compared to glass containers
- ▶ **Easy to Read Vacuum Indicator**
 - Provides visual confirmation of the vacuum inside the bottle
 - Indicates when a change of the bottle is necessary



JET-VAC 600

- Ideal for thoracentesis procedures
- Compact size is optimal for sample transfer between hospital departments



JET-VAC 1000

- Ideal for paracentesis procedures
- Excellent substitute for glass bottles



JET-VAC 1500

- Ideal for large volume paracentesis procedures
- Large size increases efficiency of procedure and reduces waste

Jet-Vac Centesis Catheters

Item No.	Catheter Size	Length	Needle Gauge	Box Qty
M7511	4F	7 cm	20G	5/Box
M7512	4F	10 cm	20G	5/Box
M7513	5F	7 cm	19G	5/Box
M7514	5F	10 cm	19G	5/Box
M7515	5F	15 cm	19G	5/Box

Jet-Vac Safety Procedure Trays (Includes Safety Needles)

Item No.	Description	Box Qty
M7552	Jet-Vac Safety Procedure Tray with 4F x 10 cm Jet-Vac Centesis Catheter	5/Box
M7553	Jet-Vac Safety Procedure Tray with 5F x 7 cm Jet-Vac Centesis Catheter	5/Box
M7554	Jet-Vac Safety Procedure Tray with 5F x 10 cm Jet-Vac Centesis Catheter	5/Box

Jet-Vac Procedure Trays

Item No.	Description	Box Qty
M7852	Jet-Vac Procedure Tray with 4F x 10 cm Jet-Vac Centesis Catheter	5/Box
M7853	Jet-Vac Procedure Tray with 5F x 7 cm Jet-Vac Centesis Catheter	5/Box
M7854	Jet-Vac Procedure Tray with 5F x 10 cm Jet-Vac Centesis Catheter	5/Box

Jet-Vac Bottles

Item No.	Item	Description	Case Qty
M7090	Jet-Vac 1000	1,000 mL Evacuated Drainage Bottle in Valve-Only Configuration	10/Case
M7092	Jet-Vac 600	600 mL Evacuated Drainage Bottle in Valve-Only Configuration	10/Case
M7093	Jet-Vac 1500	1,500 mL Evacuated Drainage Bottle in Valve-Only Configuration	10/Case
M7094	Jet-Vac 1000D	1,000 mL Evacuated Drainage Bottle with Drain Line	10/Case
M7095	Jet-Vac 1500D	1,500 mL Evacuated Drainage Bottle with Drain Line	10/Case
M7096	Jet-Vac 600D	600 mL Evacuated Drainage Bottle with Drain Line	10/Case

Accessories

Item No.	Description	Case Qty
M7005	Drainage Line Set (Includes Drainage Line with End Caps, and 5-in-1 Adapter)	10/Box
M7097	Jet-Vac Drainage Line Extension (two female luer connectors)	10/Box
M7098	Jet-Vac Large Bore Drainage Line	5/Box
M7099	Jet-Vac Drainage Line	10/Box
M7555	2000cc Drainage Bag	10/Box

U.S. Design Patent No. D818119; Canadian Design Patent No. 172153

www.pfmmedicalusa.com

pfmmedical

pfm medical, inc.
 1916 Palomar Oaks Way, Suite 150 • Carlsbad, CA 92008
 Ph: 760.758.8749 • Fx: 760.758.1167
 customerservice@pfmmedicalusa.com • www.pfmmedicalusa.com
 ©2018 PFM. All rights reserved.
 PFM Medical, the PFM logo, the color GREEN of the three-dimensional configuration of the accordion vacuum indicator and Jet-Vac are registered trademarks of pfm medical, inc.
 Chloraprep® is a registered trademark of Carefusion 2200, Inc.