


primeMIDLINE®

Power Injectable, Low-Profile Midline Catheters for reliable peripheral access of selected intravenous therapies, blood sampling and power injection of contrast media.

www.pfmmedicalusa.com

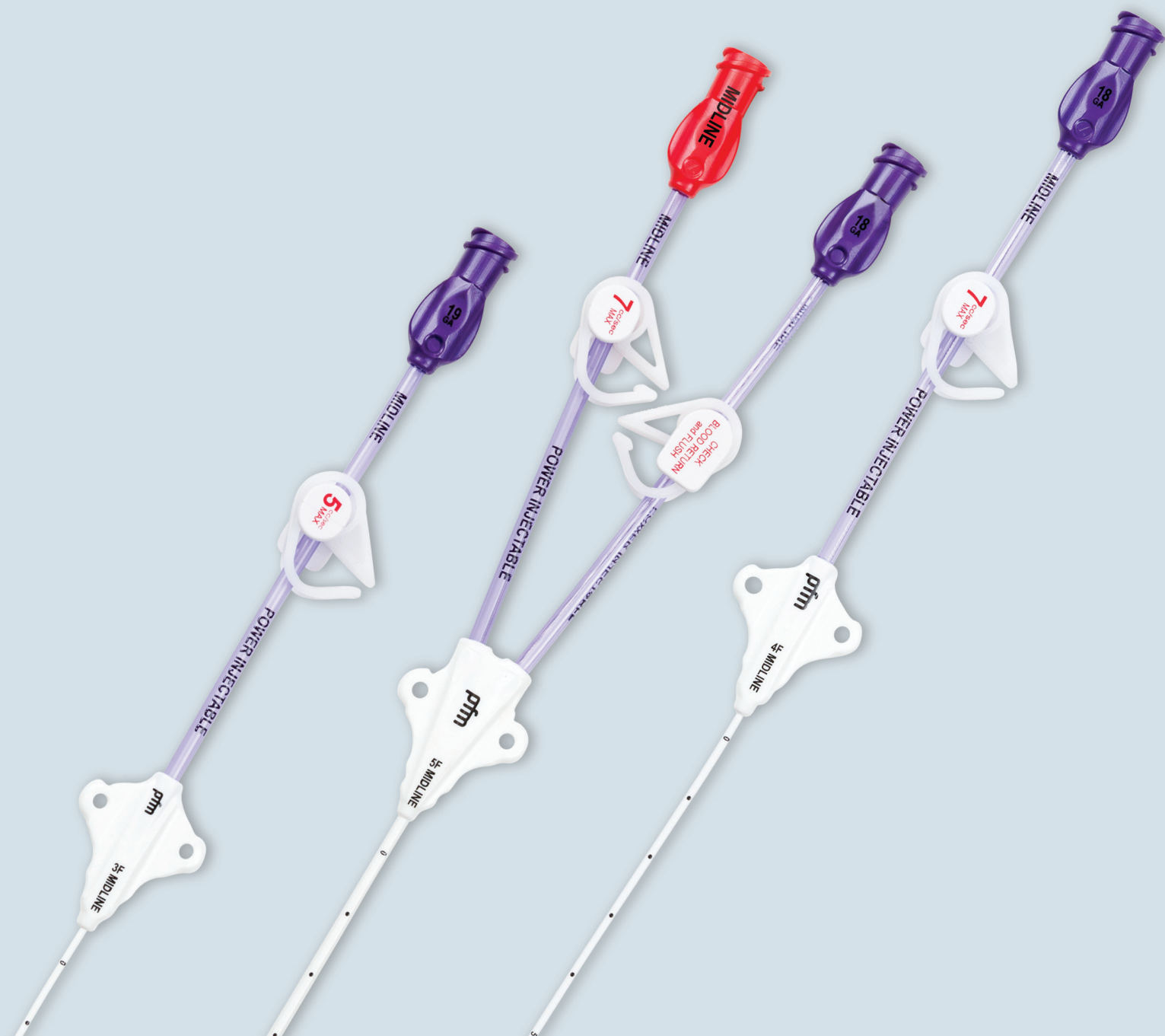

primeMIDLINE® Power Injectable, Low-Profile, Short-Term Catheters

The power injectable *primeMIDLINE* Catheter is designed to provide reliable peripheral access to the venous system. The *primeMIDLINE* Catheter combines soft polyurethane material and a low-profile luer hub for patient comfort. All lumens are power injectable to eliminate the need to select the appropriate lumen and minimize interruption of work flow. The *primeMIDLINE* Catheter offers the only high flow power injectable 3F midline catheter rated 5ml/sec.

Features & Benefits

- **Soft Polyurethane Material:**
 - Provides ease of insertion
 - Body softness to minimize risk of damage to the inner vessel wall
- **High Injection Rates Provide Optimal Imaging for Contrast Studies:**
 - 5ml/sec for 3F Single Lumen Catheter
 - 7ml/sec for 4F and 5F Single and 5F Dual Lumen catheters
- **Low-Profile Luer Hubs:**
 - Lays flat for patient comfort under occlusive dressings
 - Allows for easier grasping and twisting when changing connections
 - Resists cracking and deformation
- **3F High Performance Offering:**
 - Accommodates smaller vessels
 - Achieve high flow performance at 5ml/sec
 - Minimize thrombotic risks
- **True to Size Catheter with No Taper at the Distal End**

Kit Components

- Midline Catheter
- Tuohy Borst Side Port with Pre-Loaded Stylet
- Safety Introducer Needle 21G x 7cm with Echo Tip
- Guidewire 45cm
- Tearaway Introducer
- Syringe 10cc
- Safety Scalpel
- Tape Measure
- Neutral Needleless Valve
- Statlock® Plus Securement Device

*Statlock is a registered trademark of C.R. Bard, Inc or an affiliate

Ordering Information

REF Number	Lumens	French Size	Catheter Length (cm)	Lumen Gauge	Max CT Rate (ml/sec)	Max Pressure (psi)
PFM3SML20	Single	3F	20	19G	5.0	300
PFM4SML20	Single	4F	20	18G	7.0	300
PFM5SML20	Single	5F	20	17G	7.0	300
PFM5DML20	Dual	5F	20	18G/18G	7.0	300

www.pfmmedicalusa.com

pfmmedical

pfm medical, inc.
1916 Palomar Oaks Way, Suite 150 • Carlsbad, CA 92008
Ph: 760.758.8749 • Fx: 760.758.1167
customerservice@pfmmedicalusa.com • www.pfmmedicalusa.com

©2018 PFM. All rights reserved.
Unless otherwise noted, PFM Medical, the PFM logo, and *primeMIDLINE* are registered trademarks of pfm medical, inc.